

# EQUIPMENTS

## RIGID ENDOSCOPES

### BRV

- Custom made " Sapphire objectives" guarantee lasting
- Autoclavability, scratch resistance and large field of view
- State of art' Rod Lens systems' for superior image quality
- High tech' Tubing Material' more durable and scratch resistant and rust proof
- Optimized fiber optics provide more light



#### LAPAROSCOPES FULL HD

Description	Part Number	Dia	Length	View Angle
Laparoscope	TLPFH-05-310- 0°	5mm	310	0°
Laparoscope	TLPFH-05-310-30°	5mm	310	30°
Laparoscope	TLPFH-10-340- 0°	10mm	340	0°
Laparoscope	TLPFH-10-340-30°	10mm	340	30°

#### LAPAROSCOPES HD

Description	Part Number	Dia	Length	View Angle
Laparoscope	TLPHD-10-315- 0°	10mm	315	0°
Laparoscope	TLPHD-10-315-30°	10mm	315	30°
Laparoscope Bariatric	TLPHD-10-420- 0°	10mm	420	0°
Laparoscope Bariatric	TLPHD-10-420-30°	10mm	420	30°
Laparoscope Bariatric	TLPHD-10-420-45°	10mm	420	45°

#### LAPAROSCOPES STD

Description	Part Number	Dia	Length	View Angle
Laparoscope	TLPST-05-300- 0°	05mm	300	0°
Laparoscope	TLPST-05-300-30°	05mm	300	30°
Laparoscope	TLPST-65-300- 0°	6.5mm	300	0°
Laparoscope	TLPST-65-300-30°	6.5mm	300	30°
Laparoscope	TLPST-10-315-0°	10mm	315	0°
Laparoscope	TLPST-10-315-30°	10mm	315	30°

#### SINOSCOPE / ARTHROSCOPE STD

Description	Part Number	Dia	Length	View Angle
Sinoscope	TSINST-04-175- 0°	04mm	175	0°
Sinoscope	TSINST-04-175-30°	04mm	175	30°
Sinoscope	TSINST-04-175-45°	04mm	175	45°
Sinoscope	TSINST-04-175-70°	04mm	175	70°

#### LAPAROSCOPE / CYSTOSCOPE / HYSTEROSCOPE HD/ STD

Description	Part Number	Dia	Length	View Angle
Laparo/Cysto/Hystro	TLCHST-29-300- 0°	2.9mm	300	0°
Laparo/Cysto/Hystro	TLCHST-29-300- 30°	2.9mm	300	30°
Laparo/Cysto/Hystro	TLCHHD-29-300- 0°	2.9mm	300	0°
Laparo/Cysto/Hystro	TLCHHD-29-300-30°	2.9mm	300	30°

#### CYSTOSCOPE / HYSTEROSCOPE HD/ STD

Description	Part Number	Dia	Length	View Angle
Cystoscope	TCYSHD-04-302-0°	4mm	302	0°
Cystoscope	TCYSHD-04-302-30°	4mm	302	30°
Cystoscope	TCYSST-04-302-0°	4mm	302	0°
Cystoscope	TCYSST-04-302-30°	4mm	302	30°

#### LAPAROSCOPE SINGLE PUNCTURE [ offset ]

Description	Part Number	Dia	Length	View Angle
Laparoscope single puncture [6mm Dia Instrument channel]	TLSP-10-270-0°	10mm	270	0°

\*Telescope Box  
TTB-01/TB-02/TB-03



\*Protection Sheath  
TTPS-10-31/TPS-10-42/TPS-4-18/TPS-45-31




TTB-01 : Telescope box for Arthroscope / Sinoscopes  
TTB-02 : Telescope box for Laparoscope  
TTB-03 : Telescope box for Bariatric Scopes  
TTB-04 : Telescope box for Cystoscope/Hystroscope

TTPS-10-31 : Telescope Protection Sheath for 10mm Scopes  
TTPS-10-42 : Telescope Protection Sheath for Bariatric Laparoscope  
TTPS- 4-18 : Telescope Protection Sheath for 4mm Sinoscope  
TTPS-45-31 : Telescope Protection Sheath for 4/5mm (for 4mm Cystoscope & 5mm Laparoscope)

# ENDOSCOPIC CAMERA

## BRV CMOS HD Camera [HD Vision -C/RZ-S]

- 1/2.8" CMOS Mega Pixel Sensor
- 1080P full HD Image output, life like colour reproduction
- Image freeze
- Highest resolution of 1100 lines for better depth & sharpness
- Camera head with two Buttons (REC & WB)
- Multiple output offer flexibility for connecting to peripheral device
- Full HD recording in 1080P, Live video into USB disk
- Multiple output offer flexibility for connecting to peripheral Devices  
HD-SDI, DVI-D, VGA ,CVBS

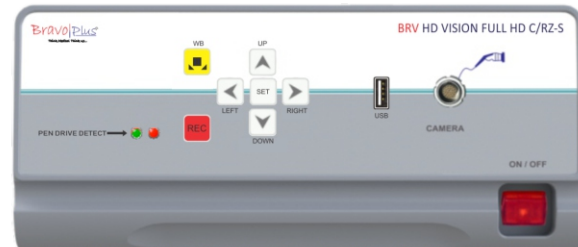

### DEVICE

Part No.	Description
BRV- C/RZ-S	BRV CMOS Full HD Camera (HD vision) with recording & all Its Standard Accessories
BRV- C- CH-S	CMOS camera head (HD Vision-CRZ)
BRV- C- CCU	CMOS control unit (CCU) (HD Vision-CRZ)
ZCO-3	Zoom coupler 20-35mm

### TECHNICAL DATA

Camera Head	1/2.8" CMOS Mega Pixel Sensor
Resolution	1920 x 1080P (1100 TV Lines )
High definition output (HD)	Full HD (16:9)
Light Sensitivity	1 Lux
Signal to Noise ratio	50 db
Shutter	1/60s-1/10000s (Automatic )
White Balance	Manual (Push)
Video output	HD-SDI, DVI-D, HDMI, CVBS
Camera Head with control button	2 Buttons
Camera Cable Length	3 meters
Power Supply	100-240 AC
Frequency	50 to 60 Hz
Maximum Power	50VA
Classification	Protection Class I, Type BF
Dimensions (WxHxD)	300 x 125 x 355 mm
Weight	4.5 Kgs.

### SUPPLIED ACCESSORIES

Part No.	Description
BNC	BNC-BNC Cable For Monitor
DVI-D	DVI-D to DVI-D Cable For Monitor
HDMI-HDMI	HDMI to HDMI cable
MC-03	Mains cord

Monitor Cable  
BNC - BNC



Monitor Cable  
DVI-D to DVI-D cable



Mains cord  
TMC-03



\*USB Pen Drive 64 GB or  
Hard Disk FAT -32



Monitor Cable  
HDMI to HDMI cable



\*Video Coupler Zoom [BRV],  
Focal length : 20-35mm  
TZCO-3



\*Endoscope Coupler  
Fixed Focus f=25mm  
TVCO-25



\*Monitor Cable  
HDMI-DVI-D cable





# HD- RECORDER

## RECORDER

- Soakable camera head
- Sony CCD high definition Chip
- High definition resolution via DVI output
- Capture the still image in SD card
- High depth of Focus
- Anti moire filter for flexible endoscope
- Window adjustment and digital zoom
- Multiple outputs for monitor compatibility



## DEVICE

Part No.	Description
BRV -169SR	BRENDO Camera 1CCD -HD [Still Recoding] with all Its Standard Accessories
BRV -169SR-CH	1CCD HD camera head (169SR)
BRV -169SR-CCU	1CCD HD unit (CCU) (169 SR)
BRV - 25	Endoscope Coupler Fixed Focus f=25mm

## SUPPLIED ACCESSORIES

Part no.	Description
BNC-BNC	Monitor Cable - BNC
DVI-D	DVI-D to DVI-D cable
SD-2	Memory Card 2gb
MC-03	Mains cord

## TECHNICAL DATA

Camera Head	1/3' CCD chip PAL
Resolution	600 Tv Lines
High definition output (HD)	1280 x 720P
Minimum illumination	0.6 lux /(F1.2)
Signal to Noise ratio	60db
Shutter	1/60s-1/10000s(automatic)
White Balance	Push to set
Window	3 level adjustable
SD card Storage	Images
Video output	Composite x2, Y/C, DVI-D, RGB
Camera Head with control button	Yes
Camera Cable Length	3 meters
Power Supply	100-240 AC
Frequency	50 Hz
Maximum Power	50VA
Classification	Protection Class I, Type BF
Dimensions (WxHxD)	300 x 65 x 355 mm
Weight	3 Kgs.

Monitor Cable - BNC  
BNC-BNC



Endoscope Coupler  
Fixed Focus f=25mm  
VCO-25



DVI-D to DVI-D cable  
DVI-D



Memory card 2GB  
SD-2



Mains cord  
MC-03



\*S-video cable  
SVHS



\* HDMI-DVID cable



\*Video Coupler Zoom [BRV],  
Focal length : 18-36mm  
ZCO-3



# MEDICAL GRADE MONITOR

## MEDICAL GRADE MONITOR

- TFT displays with wide screen format
- Back light stabilization system
- Built in multiple display mode
- True color restore
- Wide viewing angle
- Hygenic design with full glass front

### DEVICE

Part No.	Description
TLM-24	Medical Grade Lemke Monitor 24"
TLM-27	Medical Grade Lemke Monitor 27"
TMDSC-2226	Medical Grade Barco Monitor 26"
TMDSC-2326	Medical Grade Barco Monitor 26"
TU2718Q	4K Monitor 27"



Mains cord  
MC-03




TECHNICAL DATA	TLEMKE	TOMS	TBARCO
Model	LM27	U2718Q	MDSC 2226
Display size	27"	27"	26"
Support color	1.07 billion	1.07 billion	> 1 billion
Resolution	1920 x 1080(Full HD)	4K 3840 x 2160 at 60 Hz	2MP (1920 x 1080)
Pixel Pitch	0.3114mm	0.1554 mm x 0.1554 mm	0.3 mm
Brightness	700cd/m2	350 cd/m <sup>2</sup>	450 cd/m <sup>2</sup>
Contrast	1000:1	1300 to 1	1400:1
	178°	178°	178°
Response time	5ms.	5 ms.	5ms.
Input signal interface	DVI-D:24pin/3G-SDI: BNCYPbPr: BNC*3/RGBS:BNC*4 CVBS:BNC/ S-VIDEO: DIN4, VGA: D-Sub 15pin/ Fiber: FC (Optional)	HDMI, Display Port, Mini Display Port	DVI-I Single link (with HDMI support)/RGBS / YPbPr (4xBNC)/S-video (4-pin Mini DIN)/Composite video (1xBNC)/3G-SDI (1xBNC)/VGA (on DVI-I)/DisplayPort (VESA std1.1a)
Output signal interface	DVI-D:24pin/3G-SDI:BNC, CVBS:BNC, S-VIDEO:DIN4, YPbPr:BNC*3, RGBS:BNC*4		DVI-D x1 S-video (4-pin Mini DIN) Composite video (1xBNC) 3G-SDI (1xBNC)
Size	645.25(W)x392(H)x50(D)MM	611.3(W)x200.3(D)x410.2(H)MM	633(W)x411(H)x90(D)MM
Weight	9.0kgs.		12.5 kg
Weight with stand	12.5kg	9.4kg	
	VESA 100x100mm		VESA 100 mm
Power supply	DC 24V	30W	100-250 V
	55W max	99W max	120W max
Power Saving	EPA, VESA		EPA, VESA
	>50,000 hours		
	Professional Protective glass		Protective, Non-Reflective PMMA

# LED LIGHT SOURCE

## BRV LED LIGHT (BRV-180D)

- Long life > 10000 hrs and more
- Aspheric Glass provide high flux output
- Instant on/off (no warm up or cool down)
- No change in color temperature over entire color spectrum
- Digital Dimming : Intensity Control
- Does not emit ultraviolet and infrared light [3.5m W/lm]
- Reduced heat transmission to the Endoscope
- LCD Display shows brightness
- Consumes very less power
- No annual bulb change is necessary
- Comply with EMC requirements
- CCT = 5700  $\pm$  500K
- CRI = 90 [Ra>90]

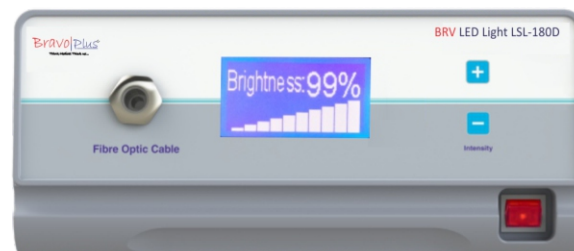

## DEVICE

Part No.	Description
BRV-180D	BRV - LED Light L180D for Laparoscopy with Light cable socket, Storz

## SUPPLIED ACCESSORIES

Part no.	Description
MC-03	Mains cord

## TECHNICAL DATA

Lamp	LED BRV-180D
Color Temperature	5700°k $\pm$ 500°k
Light guide connection	Karl Storz type, others on request
CRI	Ra>90
Lumen Maintenance	90% After 10000hrs
Ambient temperature	+10°C to +40°C in operating mode
Cooling Method	Forced Air cooling side ventilation
Frequency	50 Hz
Classification	Protection Class I, Type BF
Dimensions (W x H x D)	300 x 130 x 250 mm
Weight	6 Kgs.

Mains cord  
MC-03



\*Fibre optic cable Armoured  
FOC-523-SS



\*Fibre optic cable Non-armoured  
FXN-523-SS



## Benefit Of Aspheric Glass.

1. Aspheric plus lenses reduce magnification for a more natural appearance,
2. The decreased center thickness uses less material, resulting in less weight.
3. Traditional lenses in high-plus powers cause the wearer to look directly in the lens center for good vision



# FIBER OPTIC CABLES

## FIBRE OPTIC CABLE [FOC-523-SS]

- Projector end is heat resistance and has high duration of life
- Autoclavable at 134°C
- Fibre optic cable provides high level of light transmission
- Screwing end fittings enables connection to most light sources and instruments
- Special technology fitted ends for Xenon Safe cables. High temperature resistant
- Armoured and non armoured cables available



### DEVICE

Part No.	Description
FOC-523-SS	FOC-(230 cm, Armoured ,dia 5mm (Autoclavable, Multi-Purpose)
FXL-523-SS	FOC-(230cm, Armoured ,dia 5mm for <b>Xenon, LED &amp; Halogen Light Source</b>
FXN-523-SS	FOC-(230 cm, Light Weight, dia 5mm (Autoclavable, Multi -Purpose) for Xenon, LED & Halaogen Light Source

### TECHNICAL DATA

Length	230cms
Diameter(Active)	4.8 mm.
Diameter of Core Fibre	70 micron
Number of Fibers	3750
Compatibility Light Source End	Storz type, others on request
Compatability Scope End	Storz type, others on request
Autoclavable	134°C
Sterilizable	Soakable, ETO, Autoclavable
Classification	Protection Class I, Type BF

## Solution To Your All Needs

### Combination example

	Light Source	Light Cable	Endoscope
BRAVO PLUS			
	ALC-OM		ALE-S
Storz			
	ALC-S		ALE-A
Olympus			
	ALC-O		ALE-W
Wolf			
	ALC-W		ALE-O

# CO<sub>2</sub> INSUFFLATOR

## BRV- CO<sub>2</sub> [ AF-25L]

- High safety
- Electronic controlled gas refilling technology
- Allows easy setup and monitoring of all parameters during insufflation process
- 0.1 to 25 l/min Flow
- Micro controller based
- Inbuilt heating system
- Pressure relief system



## DEVICE

Part No.	Description
AF-25L	BRV - Insufflator for Laparoscopy, Co2 Gas-25 lpm, LCD Display With Hose Pipe

## SUPPLIED ACCESSORIES

Part No.	Description
ST-LL 3	Tube set, Insufflation reusable with luer lock
S-28/ S-16	Spanner
MC-03	Mains cord
HPT-04G	High Pressure Hose Pipe German Connection With Guage (Length = 4 Feet)

## TECHNICAL DATA

Insufflation medium	Medical Grade Carbon dioxide gas (Co2)
Display	LCD
Gas flow Maximum	25lpm
Pressure range	5-30mmHg
Power supply	100-240 V ac
Frequency	50-60Hz
Maximum Power	100VA
Classification	Protection Class 1, type Bf
Dimension (W x H x D)	300 x 138 x 320 mm
Weight	12kg

Tube set, Insufflation, reusable with luer lock  
ST-LL 3



Spanner  
S-28/ S-16



Mains cord  
MC-03



High pressure Hosepipe German Connection With Guage (Length = 4 Feet)  
HPH-04G



\* High Pressure Tubing for central gas supply- HPT-CS3



\* Filter single use only  
CF-01



\* High Pressure Hose Pipe Index Pipe Pin Index Connection With Guage (Length = 4Feet)  
HPH-PI4G



Also Available in 35L

# IRRIGATION PUMP

## BRV PUMP [IPex-1200L]


- Intelligent leak compensation for optimal view and fluid balance
- Non contact pressure measurement
- Very easy to operate
- Permanent patient monitoring
- Quick start-up time.
- Easy to load tube set
- Precise control and regulation of pressure
- Small learning curve for surgeon and nurse
- LCD Display

### DEVICE

Part No.	Description
IPex -1000	BRV - Pump [Roller Ball] With LCD Display

### SUPPLIED ACCESSORIES

Part no.	Description
ITS-03	Irrigation Tube Set with IV Needle & Luer Lock
MC-03	Mains cord

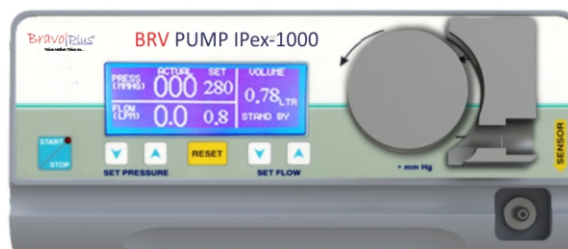

Irrigation tubing set

Mains cord

ITS-03

MC-03



### TECHNICAL DATA

Insufflation medium	Saline / Glycine
Pressure Range	5 to 400mmhg
Flow range	0.1 to 1LPM
Power Supply	170 to 240
Frequency	50 Hz
Maximum Power	150W
Classification	Protection Class I, Type BF
Dimensions (WxHxD)	300 x 320 x 138 mm
Weight	12 kgs.



### Features

- A minimally invasive technique for delivering and extracting fluid from the operating field
- Flow regulation in both tracks is carried out by the operator using a suction and flushing instrument
- Aesthetic and modern design

### Advantages of Intelligent Leak Compensation

The intelligent leak compensation automatically adjusts the inspiratory and expiratory trigger sensitivity to potential leaks.



# BRAVO SEALER

## BRV - SMART PLUS [AS-400D]

- Automatic output adjustment dependent on tissue impedance.
- Isolated RF outputs of both Monopolar & Bipolar.
- The system memory offers 10 user defined programs, plus a preset SUPPLIED program.
- Latest Micro-controller based technology used which is upgrade in features.
- Two Bipolar outputs.
- Memorization of the last operating condition.
- Portable and user friendly design.
- Lesser charring and drag.
- Maximum output is 400 cutting.
- Endocut function for polyp cutting

DEVICE	Description
BRV-400D	BRV - Smart Plus Digital [Electrosurgical Unit] With All Its Standard Accessories 400 W ( Monopolar & Bipolar)

### Application

- Gastrectomy
- Appendectomy
- Colectomy
- Gastric Bypass
- Vaginal Hysterctomy
- Colon Resection
- Redical Prostatectomy
- Whipple
- Salpingoectomy
- Lap- Cholecystomy
- TLH
- Abdominal Hysterctomy
- Nephretomy
- Liver Resection
- Splenectomy



### Monopolar Cut Modes

Types of Current	Max-Power Watt	Rated Load Ohm	Crest Factor
Pure Cut	400W	500 Ohms	CF 1.5
Blend Cut	250W	500 Ohms	CF 1.5 - 1.8
Endo Cut	100W	500 Ohms	CF Variable

### Bipolar Modes

Types of Current	Max- Power Watt	Rated Load Ohm	Crest Factor
	70W	100 Ohms	CF 1.5
	90W	100 Ohms	CF 1.5

### Monopolar Coag Modes

Types of Current	Max-Power Watt	Rated Load Ohm	Crest Factor
Spray	120W	500 Ohms	CF 7.0
Fulguration	150W	500 Ohms	CF 5.5
Dessicate	120W	500 Ohms	CF 4.5

### Safe Seal Modes

Lap	4Level
Open	4Level

SS Patient plate  
SSP-02



Patient plate cable  
PPC-02



Tri Pedal Footswitch  
FSC-03-AS



Mains cord  
MC-03



Monopolar chuck handle  
MCH-01



\*Bipolar cable  
BCS-04-03-VMF



\*Open bipolar tweezer  
OBT-07



\*Silicon Patient plate  
SPP-02



Product Picture Is For Reference

www.bravoplus.in

TRYKER

# ULTRA GRASSO

## Ultra Grasso

Ultra Grasso simultaneously seal and cut vessels up to and including 7mm\*\* in size with minimal thermal spread using combination of ultrasonic and advanced bipolar energies delivered through a single multi functional instrument

Ultra Grasso designed to provides precise, controlled dissection and available in bipolar coagulation without compromising grasping ability.

### DEVICE

ABEAT-300 USC

BRV- Ultra Grasso With  
all Standard Accessories

### SUPPLIED ACCESSORIES

Part No.	Description
ABEAT-300G	Ultrasonic Generator
BRV- AT-36-C55	Ultrasonic Transducer
BRV- ALS-36HA-C55	BRV- Laparoscopic Shears (L=36cm )
BRV- FSABU-02	Dual Pedal-footswitch for Altrabeat

ABEAT-300 USBUC Ultra Grasso with  
bipolar all Standard Accessories

### SUPPLIED ACCESSORIES

Part No.	Description
BRV- 300G	BRV- Ultrasonic Generator
AA-150CB	BRV- Accuseal Plus (ESU-bipolar)
BRV- AURT-36-C55	BRV- U&R Transducer
BRV - ALBS-36HA-C55	BRV - Laparoscopic Bi- Shears (U&R) (L=36cm)
BRV- FSABUR-03	Tri Pedal- footswitch for BRV U&R Lap Bi-Clamp Grasper
NLPG-05-30-31-SA	Lap Bi-Clamp Cutter
NLCP-05-30-31-SA	Bipolar Cable x 2
BCS-04-03-VMF	Tri Pedal-Footswitch for Accuseal (ESU)
FSC-03-AAP	

### Ultrasonic Shears



### Ultrabeat Bi-Shears (U&R)




### TECHNICAL DATA

Operating Temperature	10-40°C
Driver Frequency	55.5kHz
Scalpel Vibration	50-120um
In Put Power	AC110/220V 50Hz
Maximum Output	80W
Classification	Protection class I, Type BF
Dimensions (WxHxD)	565 x 245 x 565mm
Weight	10kg

### Application

Ultrasonic scalpel has the advantages of less bleeding, less damage to surrounding tissue, early healing, etc. The applications of ultrasound in medical surgery include incision and hemostasis, phacoemulsification, tumor aspiration, fat aspiration.


BRV Laparoscopic  
Shear  
BRV- ALS-36HA-C55#

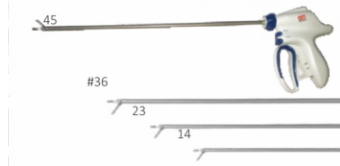

BRV Transducer  
BRV- AT-36-C55#



BRV Laparoscopic  
Bi- Shears (U&R)(L=36cm  
BRV-ALBS-36HA-C55

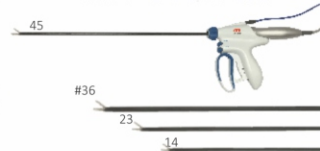

BRV (U&R) Transducer  
BRV-AURT-36-C55



Dual-foot switch for  
BRV - U&R  
BRV-FSABUR-02



Tri Pedal- foot switch for  
BRV - U&R  
BRV-FSABUR-03



\*Ultrasonic Handpiece  
(focus)  
BRV-F-9B



\*Ultrasonic Transducer  
BRV-AT-F-9B



Lap Bi-Champ Grasper/ Cutter  
NLPG-05-30-31-SA/  
NLCP-05-30-31-SA



# BRENDOTRAINER

## BRENDOTRAINER

- Basic hand eye co-ordination and to learn minimal access surgery techniques (basic and advance)
- To learn intra and extra corporeal suturing and knot tying
- Video Scopic Endo trainers are the most user friendly trainer which can be attached to any Monitor/TV.
- Built illumination & high resolution camera
- 5 Port with silicon rubber gasket.

## DEVICE

Part No.	Description
BRV-1-C	BRV Trainer (Dome Shaped) with Camera, Light & Its Complete Accessories
BRV-1-LK	BRV Trainer (Dome Shape) with Camera & Light Only
BRV-1-L	BRV Trainer (Dome Shape) with Light
BRV-1	BRV Trainer (Flat Shape) ( Box Only)
BRV	BRV Trainer Accessories Set



## TECHNICAL DATA

Power supply	100-240 V ac
Frequency	50-60Hz
Maximum Power	50W
Classification	Protection Class 1, type Bf
Dimension (W x H X D)	450mmX190mmX295mm
Weight	12kg

## SUPPLIED ACCESSORIES - BRV- 1-C

Part no.	Description
BRV -05-30-M3C	Plain Grasper 5mm –NDM
BRV -05-30-M3C	Dissector Maryland 5mm - NDM
BRV -05-30-M3C	Curve Scissor 5mm - NDM
BRV -05-30	Needle Holder T.C. Axial Type straight 5mm - NDM
BRV -1	Mains 2 Pin Cable (for DSPT-1, DSPT-1-C)
RCA to RCA	Video cable
BRV -01-40-12	knot Pusher (Ø=1.2mm) with Reducer
BRV -4X4	Silicon Suture Pad

Bhandarkar's knot Pusher  
BRV -05-40-12



Plain Grasper 5mm  
BRV-05-30-M3C



Dissector Maryland 5mm  
BRV-05-30-M3C



Curved Scissor 5mm  
BRV-05-30-M3C



Needle Holder TC Axial Straight/  
Curve 5mm\*\*  
BRV-05-30/ BRV-05-30



Mains 2 Pins Cable  
BRV-1



Video cable  
RCA to RCA



Silicon Suture Pad  
BRV-4X4



\* Washer  
WET-05/WET - 10



One Of The Best Product To Train Yourself



# TROLLEY

## Specifications:


- Overall Approx Dimensions : 950 mm x 590 x 1600 mm
- Top Shelf Size : 720mm x 440mm
- Other 4 Shelves : 575mm x 440mm
- Drawer at bottom : 575mm x 590mm
- Five M.S. Shelves
- All Shelves Are Height Adjustable
- 4 x 125 mm high grade synthetic castors
- Brake mechanism ( on two castors)
- Stainless steel IV rod and power coated gas cylinder cage
- Two stainless steel push handles (Grade : 304)
- Suitable for keeping monitor, light source, camera unit, VCR, Printer, stabilizer etc.
- 1 MCB Installed for safety
- 5 x High Grade Switches
- CNC Punched And CNC Bent
- Pretreated in 7 tank process and Epoxy Power Coated
- Anti-static wheels

**EASY TO  
INSTALL &  
PROTECT  
EQUIPMENT**



**One Of The Best Trolley Designed For Laparoscopy**

# NASOFIBROSCOPE

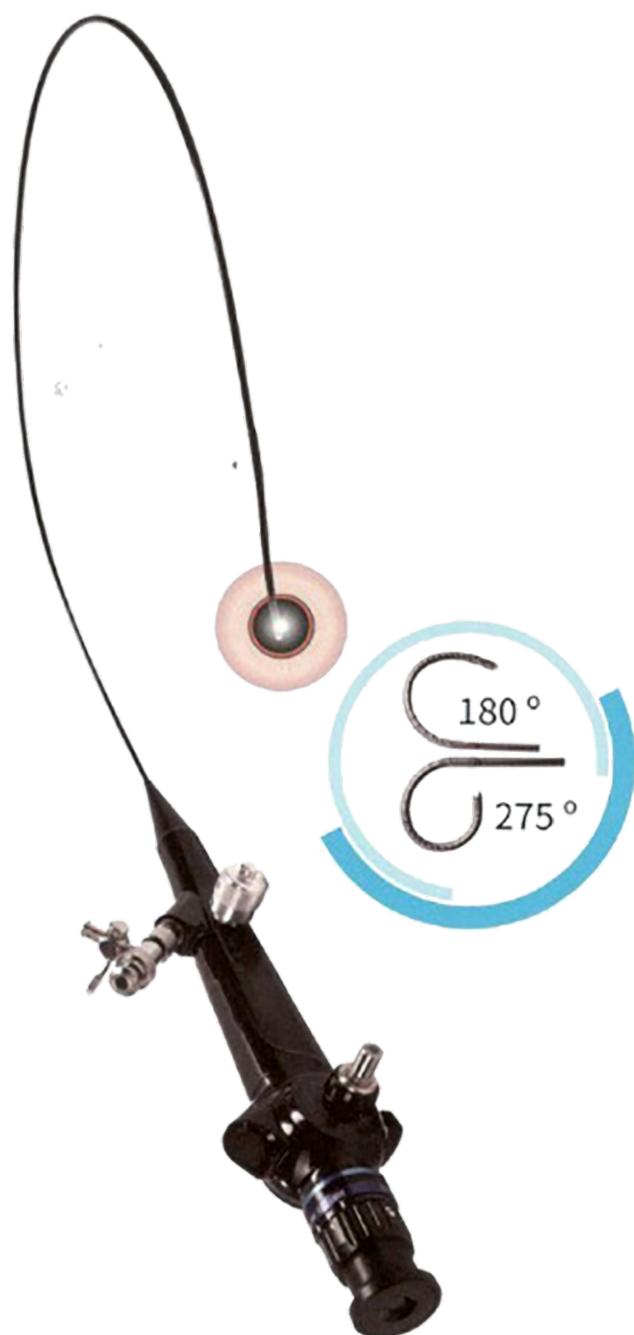

- Bright observation image for easy diagnosis
- 2.8mm extra-slim outer diameter, waterproof, easy and smooth insertion
- 1.2mm instrument channel, compatible for all kinds of guidewire passing by
- Flexible up and down angulation at 180° and 275° respectively, easily satisfies urologists examination, diagnosis and treatment requirements in the urinary tract and upper & lower renal calyx
- 40x eyepiece, imported high-definition image bundle and light bundle
- Working together with high brightness lamp light source, video system and endoscopy workstation, not only greatly facilitating operation, but also enabling typical image to be archived and providing reference for future comparison
- Suitable for all kinds of video system with flexible or rigid scope lens adapters

field of view	Depth of view Field	Range of Angulation	Outer Diameter of Distal End	Outer Diameter of Insertion Tube	Inner Diameter of Instrument Channel	Working Length	Waterproof
Forward Viewing	3mm-50mm Fixed Focus	180° 275°	2.8mm	2.8mm	1.2mm	730mm	

Specification design and accessories are subject to change without any notice on obligation on the part of manufacture